

OR*SPECIFIC*

MEGA

instrument table



Maximize Space

Upper shelf creates additional space without compromising visibility



Visibility

Easily arrange & organize instrument trays without the need for stacking



Time

Improve setup efficiency & workflow during large surgical cases



Drape

Create a sterile field with our patented one-piece drape.

Standard or heavy-duty


Protected by U.S. Patent No. 6,019,102





adjustable **BIG CASE Back Tables**

OR Specific's adjustable back tables accommodate variable staff heights ranging from 4'11" to 6'6". The main surface of the table may be adjusted intraoperatively from a minimum height of 29" to a maximum height of 38".

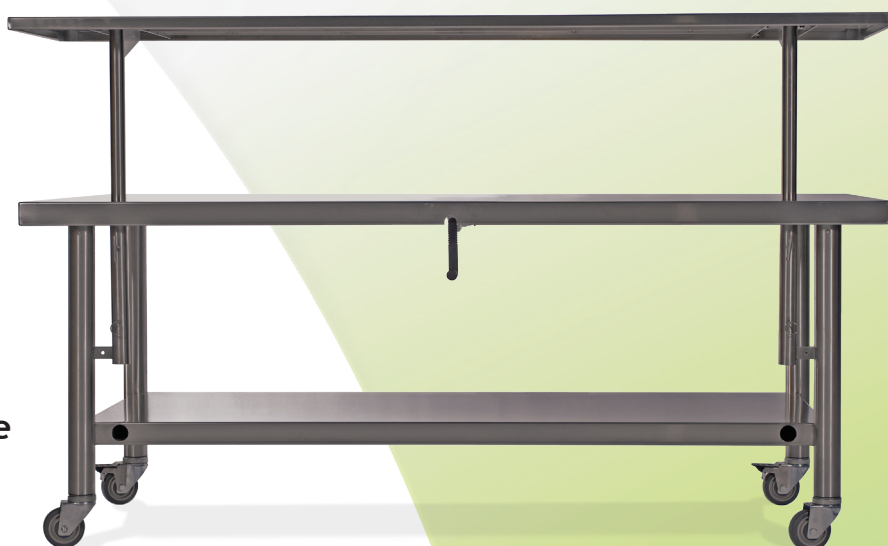

428-ADJ Big Case 5' Adjustable

Table Dimensions: 60"L x 30"D x 29"H

Adjustable: 29" LOW to 38" HIGH


Shelf Dimensions: 60"L x 21½"D




Two locking casters on front legs of table hold table stationary.



Table shown with upper tier in storage position

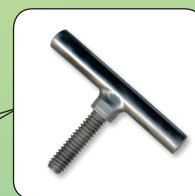

Secure the upper tier with T-bolts once raised and positioned.

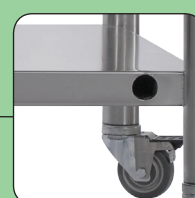

Post hole for storing upper tier when not in use.

MEGA

instrument table

ADJUSTABLE TABLE SPECIFICATIONS

Adjustable Big Case Back Table models may be raised or lowered intraoperatively through the use of a hand crank.

- **Adjustable Main Surface**
29" minimum height
38" maximum height
- **Upper Shelf Set Positions**
12" above main surface
15" above main surface
- **Weight Capacity**
150 lbs. upper shelf
200 lbs. main surface
200 lbs. storage shelf
— weight uniformly distributed
- **18 gauge stainless steel**
- **Four 3½" casters**
- **Two front locking casters**



429-ADJ Big Case 6' Adjustable


Table Dimensions: 72"L x 30"D x 29"H

Adjustable: 29" LOW to 38" HIGH


Shelf Dimensions: 72"L x 21½"D



Raise or lower the table by using the hand crank at back of table.

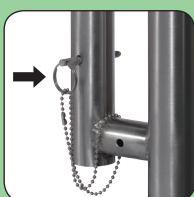

Safety locking pin used for securing upper tier shelf.

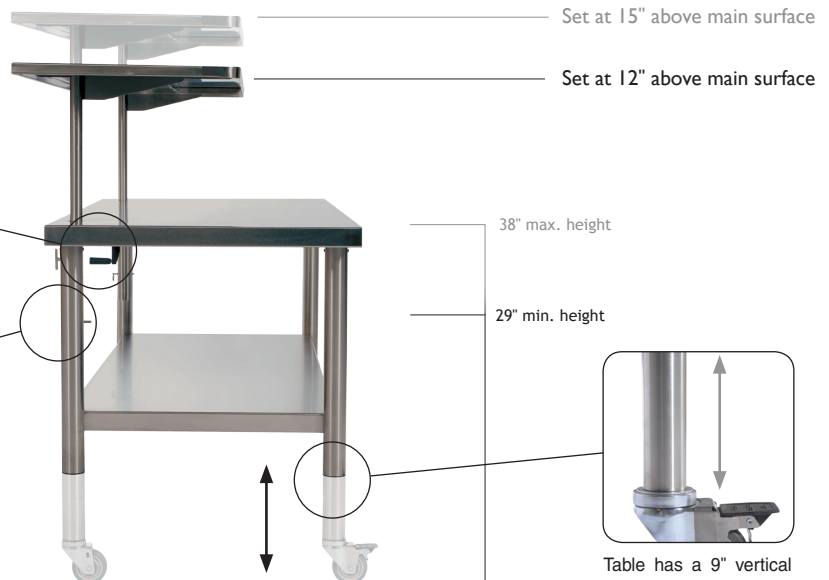



MEGA

instrument table

Two-Tier Table Drapes (Protected by U.S. Patent No. 6,019,102)

- Unique single piece Heavy Duty disposable drapes:
- Three layer system including extra absorbent center layer - 28 mm thick.
- One piece drape easily covers both upper and main table surfaces.
- The clear plastic backdrop protects the sterile field while allowing light through to the main table surface.

*Guaranteed True-Fit Drapes
for OR Specific's Big Case Back Tables*




Heavy-Duty Drape

Part No.	Description	Dimensions	Drap Part No.	Draps Per Case
428-ADJ	5' Back Table with removable 21½"D second tier self	60"L x 30"D x 29"H (height adjustment from 29" to 38")	420-HD-S	10
429-ADJ	6' Back Table with removable 21½"D second tier self	72"L x 30"D x 29"H (height adjustment from 29" to 38")	419-HD-S	10



Draping Video DVD Demo

Request a copy of the features, benefits and draping instructions for the Big Case Back Table. CD-ROM or DVD. Also available as a download at www.orspecific.com.

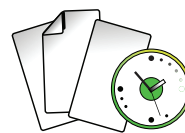

Efficiency Report

Request a copy of the time saving benefits for the Big Case Back Table from your local rep. Also available as a PDF download from our website at www.orspecific.com.



TEL +1 360.567.1972

FAX +1360.696.1700

4000 SE Columbia Way
Vancouver, WA 98661 USA

www.orspecific.com



Distributed by



Cutting Edge Innovation

3/17 Rodborough Rd, Frenchs Forest NSW 2086
PO Box 6052, Frenchs Forest DC NSW 1640, Australia
Customer Service T: 1300 665 884 F: 1300 665 886
www.surgicalspecialties.com.au

SS052